



# PrimaConnex® Surgical Drill Chart

## Tapered Implants • Sequence and Depth Marking Information (Sequence for a Ø4.1x13mm, RD implant shown)\*

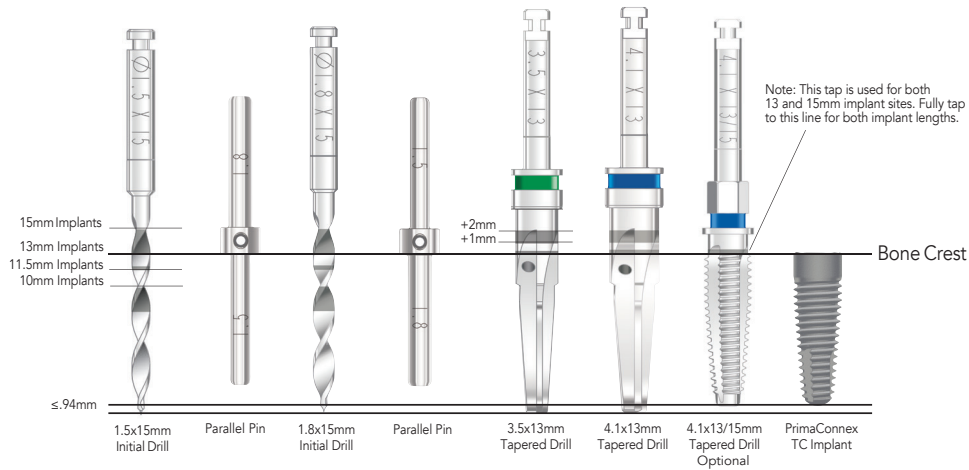

*Note: Drill bands reference the insertion depth of the corresponding implant and do not indicate actual millimeter drill length. Due to the drill design, Keystone Dental Prima Tapered Implant Drills are a maximum of .94mm longer than the insertion depth of the implant.*

## Straight Implants • Sequence and Depth Marking Information (Sequence for a Ø4.0x13mm, RD implant shown)\*

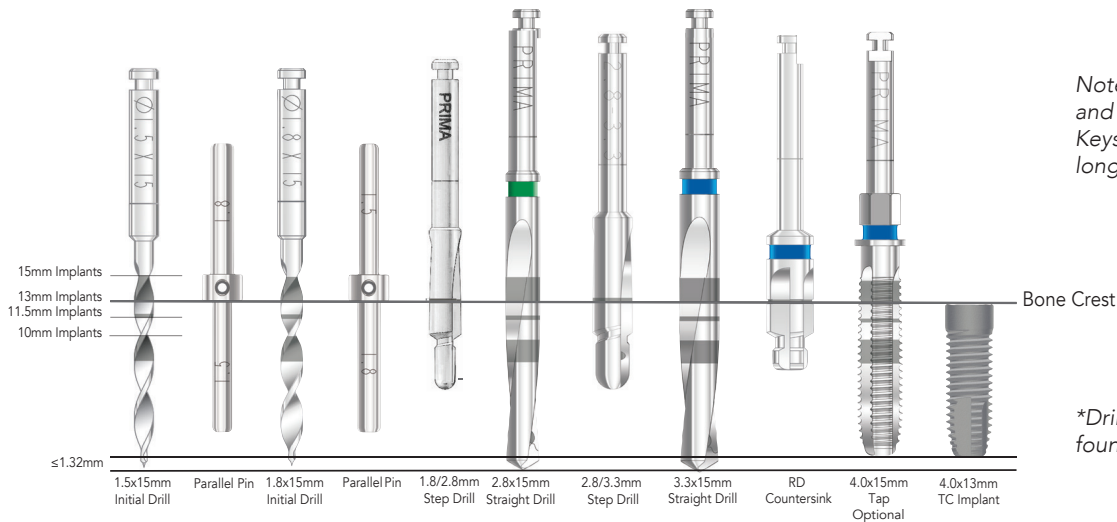

*Note: Drill bands reference the insertion depth of the corresponding implant and do not indicate actual millimeter drill length. Due to the drill design, Keystone Dental Prima Straight Implant Drills are a maximum of 1.32mm longer than the insertion depth of the implant.*

*\*Drill sequence charts for all Keystone Dental Prima Implant sizes can be found in the Keystone Dental Prima Implant System Surgical Manual.*